# **Kreatech™ FISH probes MYC (8q24) Break – XL for BOND**

**Catalog No: KBI-XL006**

- EN** - Instructions for use
- FR** - Mode d'emploi
- DE** - Gebrauchsanleitung
- IT** - Istruzioni per l'uso
- ES** - Instrucciones de uso
- PT** - Instruções de Utilização
- NL** - Gebruiksaanwijzing
- SV** - Instruktioner för användning
- DA** - Brugsanvisning
- TR** - Kullanma talimatları

## Table of Contents

<b>EN</b>	- Instructions for use	5
<b>FR</b>	- Mode d'emploi	17
<b>DE</b>	- Gebrauchsanleitung	31
<b>IT</b>	- Istruzioni per l'uso	45
<b>ES</b>	- Instrucciones de uso	59
<b>PT</b>	- Instruções de Utilização	71
<b>NL</b>	- Gebruiksaanwijzing	83
<b>SV</b>	- Instruktioner för användning	95
<b>DA</b>	- Brugsanvisning	107
<b>TR</b>	- Kullanma talimatları	119

# Kreatech™ FISH probes

## MYC (8q24) Break – XL for BOND

Catalog No: KBI-XL006

### Intended Use

For *in vitro* diagnostic use

MYC (8q24) Break – XL for BOND is intended for use on lymphoma formalin-fixed, paraffin-embedded (FFPE) tissue sections. The probes are to be used in combination with the BOND FISH Kit (DS9636) on an automated BOND system (includes BOND-MAX system and BOND-III system).

The MYC (8q24) Break – XL for BOND FISH probe detects genomic translocations involving the MYC gene in FFPE tissues. The probe has been validated on lymphoma tissues.

The clinical interpretation of any staining or its absence should be complemented by morphological studies and proper controls and should be evaluated within the context of the patient's clinical history and other diagnostic tests by a qualified pathologist.

### Summary and Explanation

Kreatech FISH probes are used to perform a fluorescence *in situ* hybridization (FISH) staining procedure on formalin-fixed, paraffin-embedded tissues. Following appropriate pretreatment, hybridization of Kreatech FISH probes and post-hybridization stringent washing using the BOND FISH kit, tissue sections are to be manually dehydrated and mounted with DAPI. Results are interpreted by fluorescence microscopy using the recommended filters with the appropriate wavelengths.

## Reagents Provided


Product Description	Product No.	Volume	Composition
MYC (8q24) Proximal – XL	pKBI-XL006G	1 ml, 10x concentrate	PlatinumBright™ 495 fluorophore-labeled DNA probes in hybridization buffer. Hybridization buffer contains < 50% Formamide
MYC (8q24) Distal – XL	pKBI-XL006R	1 ml, 10x concentrate	PlatinumBright™ 550 fluorophore-labeled DNA probes in hybridization buffer. Hybridization buffer contains < 50% Formamide

MYC (8q24) Proximal – XL is optimized to detect the genomic regions proximal to break points in the MYC gene region.

MYC (8q24) Distal – XL is optimized to detect the genomic regions distal to break points in the MYC gene region.

Both probes combined are used to detect translocations involving the MYC gene at 8q24.

**Table 1:** Chromosome map showing the size of the probes in kilobasepair (kB), chromosome band, and genomic position defined by STS markers<sup>1</sup>. The arrow adjacent to the gene indicates the transcription direction of the gene. In green MYC (8q24) Proximal – XL. In red MYC (8q24) Distal – XL. See the Leica website for more information on chromosome maps.

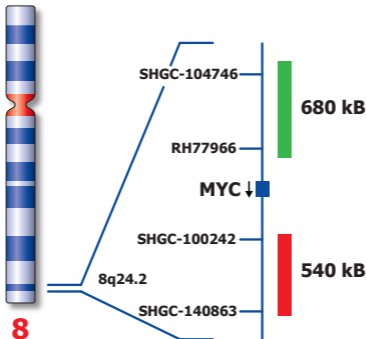- 1) Maps and locations are based on UCSC genome browser, build GRCh37/hg19, Feb.2009.

### Materials Required But Not Provided

Refer to “Using BOND Reagents” in your BOND user documentation for a complete list of materials required for specimen treatment and *in situ* hybridization using the BOND system (includes BOND-MAX system and BOND-III system).

In addition to the generic BOND reagents, the following reagents are specifically required for the analysis of suitable patient material using Kreatech FISH probes on the BOND system (includes BOND-MAX system and BOND-III system).

- BOND FISH kit (DS9636)

Kreatech FISH probes are run on a BOND system (includes BOND-MAX system and BOND-III system) in combination with the BOND FISH kit (DS9636). The BOND FISH kit in combination with \*FISH Protocol D, preinstalled on the BOND system with BDZ v62 or higher, is necessary for FISH to be run on the BOND system. This protocol is designed as a post-hybridization wash protocol, to reduce non-specific hybridization of nucleic acid probes. The BOND FISH kit contains the post-hybridization wash buffer and a single kit is sufficient for 60 tests.

- BOND Hybridization solution (AR9037)

The probes are to be diluted to the preferred concentration using the BOND Hybridization solution.

- BOND Enzyme Pretreatment kit (AR9551)

Enzyme digestion in combination with heat pretreatment is required for optimal FISH results.

- DAPI (e.g. LK-095A or LK-096A)

DAPI counterstaining is not automatically performed as part of the FISH protocol on BOND and is to be done manually using the preferred concentration of DAPI counterstain / antifade.

## Specimens

Kreatech FISH probes have been tested on FFPE tissues which were cut at ~5 µm. Tissue blocks and cut slides should be stored according to the recommended storage and handling conditions for optimal results, see Leica website.

Baking can be done either on the BOND using the \*Bake and Dewax protocol or separately in an oven or incubator. Recommended baking temperature and time is up to 16 hrs at 56 °C, or 60 - 120 min at 80 °C.

## Instructions for use: Dilution and Mixing

Kreatech FISH probes are provided in a 10x concentrated format and are to be combined and diluted to a 1x final concentration using the BOND Hybridization Solution (AR9037).

Reagent	Volume used (recommended)
Probe 1 (10x) - Red	1.0 ml
Probe 2 (10x) - Green	1.0 ml
BOND Hybridization solution (AR9037)	8.0 ml
Total volume (1x final concentration)	10 ml
# Tests	Up to 40 tests

After dilution of the probe, the mixture is to be transferred to an appropriate BOND container (titration container, 7 ml open container or 30 ml open container, depending on prepared volumes) and registered as "DNA probe" in the BOND.

When preparing dilutions for titration, ensure that allowance is made for the dead volume of the appropriate BOND Open Container (for values refer to your BOND user documentation). The exact number of tests from a vial depends on the dead volume of the used container and the calibration of the BOND system.

See Troubleshooting section for use of final concentrations other than 1x.

Refer to “Reagent Management” in your BOND user documentation for information about registering new reagents on the BOND system (includes BOND-MAX system and BOND-III system).

### Instructions for use: Pretreatment and Protocols

Kreatech FISH probes have been validated using the BOND FISH kit, BOND ancillary reagents and the BOND pretreatment and staining protocols shown in Table 2. The final protocol may vary according to the tissue used by the user. Any deviations from the prescribed protocols must be validated by the user.

**Table 2:** Protocols for running Kreatech FISH probes

BOND Protocol Section	Protocol Name	Protocol Type
Staining – ISH Detection	*FISH Protocol D	Post Hybridization Wash
Prestaining – Preparation	*Dewax	Removal of paraffin
Prestaining – Heat Pretreatment	HIER 25 min with ER2 (100) <sup>(1)</sup>	Heat pretreatment, Epitope retrieval
Prestaining – Enzyme Pretreatment	*Enzyme 5 for 20 min <sup>(1,2)</sup>	Enzymatic digestion
Prestaining – ISH denaturation	*D10	Denaturation 10 min.
Prestaining - ISH detection	*ISH Hybridization (12Hr)	Hybridization 12 hr.

- (1) Not a predefined protocol. For creating a new protocol refer to “Protocols” in your BOND user documentation.
- (2) Enzyme 5 is derived from Enzyme Concentrate 1 available in the BOND Enzyme Pretreatment Kit (AR9551). See Troubleshooting section below for specifics.

### Signal patterns

MYC (8q24) Break – XL for BOND can be used to detect translocations involving the MYC gene in FFPE tissues. The probe has been validated on lymphoma tissues using the protocols shown in table 2.

MYC (8q24) Break – XL for BOND is designed as a dual-color assay to detect translocations of the MYC gene region at 8q24. The normal signal pattern shows two red/green fusion signals (2F), while a translocation involving the MYC gene results in one fusion signal for the normal chromosome 8 and one red and one green signal for the translocation (1F1R1G).

Follow-up testing is advised if other, atypical patterns are observed.

	Normal Signal Pattern	Translocation involving MYC gene at 8q24
Expected Pattern	2F	1F1R1G

**Disclaimer:** The performance of MYC (8q24) Break – XL for BOND has been established using the procedures described in this package insert. Modifications to these procedures may alter the performance of the assay. Any deviations from the instructed protocols must be validated by the user.

### Storage and Stability

Protect probes from strong and direct light.

Store at 2 - 8 °C. The unopened product is stable under these conditions up to the expiry date indicated on the vial label.

Return undiluted reagent to 2 - 8 °C immediately after use. Diluted/ mixed probe reagent stored in a BOND container should be returned to 2 - 8 °C after being run.

After probe dilution, probes mixtures are stable for up to 30 days. After this time period the diluted mixture may result in suboptimal results.

There are no obvious signs that indicate contamination and/or instability of the probe.

Storage conditions other than those specified above must be verified by the user.

### **Troubleshooting**

Pretreatment (overdigestion): The Instructions For Use of the BOND Enzyme Pretreatment kit (AR9551) recommends a dilution of 1 drop of Enzyme Concentrate 1 into 7 ml of Diluent to be registered as \*Enzyme 1 in the BOND software.

This \*Enzyme 1 may be too strong for certain specimens when used in combination with the Heat Induced Epitope Retrieval (HIER) and results in overdigested tissue. A secondary dilution of the \*Enzyme 1 (up to 1:300 dilution) is recommended and a separate registration of the reagent on the BOND system as \*Enzyme 5 is recommended to prevent confusion. Alternatively, a direct dilution of up to 1:3000 of Enzyme concentrate 1 can be used instead of the serial dilution of \*Enzyme 1.

Probe dilution and Signal intensity: A 1x final probe concentration is recommended for use but optimal probe balance is dependent on factors such as tissue age, tissue fixation, pretreatment protocol, and personal preference. Optimal concentrations and probe ratios may vary and are to be determined empirically by the end-user.

Reagent	Volume used (recommended)	Example: Reduce Green signal strength	
Probe 1 (10x) - Red	1.0 ml	Red at 1x:	1.0 ml
Probe 2 (10x) - Green	1.0 ml	Green at 0.5x:	0.5 ml
BOND Hybridization solution (AR9037)	8.0 ml		8.5 ml
Total volume at 1x final concentration	10 ml		10 ml
Number of tests	up to 40 tests		up to 40 tests

Contact your local distributor or the customer care of regional office of Leica Biosystems to report any unusual staining.

### Product Specific Limitations

Users who deviate from recommended test procedures must accept responsibility for interpretation of patient results under these circumstances. Protocol times may vary due to variation in tissue type, fixation and processing. In addition, BOND Enzyme concentration and incubation time may require optimization depending on tissue type, processing and fixation conditions.

### Further Information

Further information on *in situ* hybridization with BOND reagents can be found in your BOND user documentation.

### Warnings and Precautions

- This product is intended for *in vitro* diagnostic use only.
- For laboratory / professional-use only.

- Kreatech FISH probes are intended for use on BOND systems.
- Do not use reagents after the expiry date given on the package label.
- Wear suitable protective clothing, gloves and eye/face protection.
- Kreatech FISH probes contain formamide and should therefore be registered in the BOND software as hazardous waste.
- Formamide may have adverse effects upon human reproduction. Avoid inhalation, ingestion or contact with unprotected skin. If contact with skin occurs, wash with copious amounts of water for at least 15 minutes.

**DANGER****FORMAMIDE****Hazard statement****H360D** May damage the unborn child.**Precautionary statements****P201** Obtain special instructions before use.**P280** Wear protective gloves/protective clothing/eye protection/face protection.**P308+P313** If exposed or concerned: Get medical advice/attention.**P405** Store locked up.**P501** Dispose of contents/container to authorised chemical landfill or if organic to high temperature incineration.

Consult the SDS for full details. To obtain a copy of the Safety Data Sheet contact your local distributor or regional office of Leica Biosystems. Safety Data Sheets are available on the Leica Biosystems website, [www.LeicaBiosystems.com](http://www.LeicaBiosystems.com)

- Specimens, before and after fixation and all materials exposed to them, as well as leftover probe volume, should be handled as if capable of transmitting infection and disposed of with proper precautions.

- Never pipette reagents by mouth and avoid contacting the skin and mucous membranes with reagents or specimens.
- If reagents or specimens come into contact with sensitive areas, wash with copious amounts of water. Seek medical advice.
- Clean up spillages immediately using suitable disposable materials.
- All materials must be disposed of in accordance with local legislation.
- Retrieval steps, incubation times or temperatures, other than those specified, may give erroneous results. Any such change must be validated by the user.

### Bibliography

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide.

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## Bibliographie

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue ; proposed guideline. 1991 ; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In : Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York : Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide (Guide de l'utilisateur de FISH pour BOND).

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### Bibliographie

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide (Bedienungsanleitung zu FISH auf BOND).

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## Bibliografia

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide (Guida per l'utente di FISH su BOND)

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1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide (Guía del usuario para FISH en el sistema BOND).

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## Bibliografia

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide (Guia do utilizador FISH no BOND)

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## Bibliografi

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide (FISH on BOND Bruksanvisning).

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## Bibliografi

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Bestillingsnr. M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. I: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, s.1 - 13.
4. FISH on BOND User Guide (Brugervejledning til FISH på BOND)

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### Kaynakça

1. Clinical Laboratory Improvement Amendments of 1988, Final Rule 57 FR 7163 February 28, 1992.
2. Villanova PA. National Committee for Clinical Laboratory Standards (NCCLS). Protection of laboratory workers from infectious diseases transmitted by blood and tissue; proposed guideline. 1991; 7(9). Order code M29-P.
3. Wilkinson DG. The theory and practice of *in situ* hybridization. In: Wilkinson DG. (ed.) *in situ* Hybridization. A practical approach. 2nd Edition. New York: Oxford University Press, 1998, pp.1 - 13.
4. FISH on BOND User Guide (FISH on BOND Kullanıcı Kılavuzu)

DANGER

FORMAMIDE

2°C

8°C  
long term  
storage

Kreatech Biotechnology B.V.  
Vlierweg 20  
1032 LG Amsterdam  
The Netherlands  
[www.LeicaBiosystems.com](http://www.LeicaBiosystems.com)

For detailed information about European sales offices,  
please visit our website: [LeicaBiosystems.com](http://LeicaBiosystems.com)

**Patents:** These products or the use of these products is subject to proprietary rights. The probes in these products are produced using the Janssen Diagnostics, LLC REPEAT-FREE™ technology and labeled with the Universal Linkage System (ULS™).

The ULS™ technology and products are covered by one or more of the following US patents 5,580,990; 5,714,327; 5,985,566; 6,133,038; US RE 40,557E, 6,797,818; 7,217,813. KREATECH is a trade name of KREATECH Biotechnology BV. ULS™, PlatinumBright™, are trademarks of KREATECH Biotechnology BV. REPEAT-FREE™ is a trademark of Janssen Diagnostics, LLC.

KBI-XL006 Rev B 12/02/2016 · Copyright © by Leica Biosystems, Amsterdam, 2014. Subject to modifications. LEICA and the Leica Logo are registered trademarks of Leica Microsystems IR GmbH.