

Kreatech™ FISH probes Product Information Sheet

18Q003V495

BCL2 (18q21) Proximal - XL


Green

1 ml

DANGER



FORMAMIDE

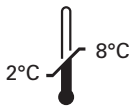


Kreatech Biotechnology B.V.
Vlietweg 20
1032 LG Amsterdam
The Netherlands
www.LeicaBiosystems.com

Analyte Specific Reagent

Analytical and performance characteristics are not established

PI-18Q003V495_D1.2

Published October 2016

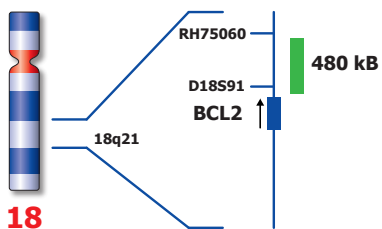

Not to scale

18Q003V495

Kreatech™ BCL2 (18q21) Proximal -XL FISH probe

Intended use: Intended for use in a diagnostic application for identification and quantification of an individual chemical substance or ligand in biological specimens.

Introduction: The **BCL2 (18q21) Proximal - XL** FISH probe is optimized to detect the genomic regions proximal to the BCL2 breakpoint(s) at 18q21.

Critical region (green): The **BCL2 (18q21) Proximal - XL** FISH probe is direct-labeled with PlatinumBright™495.

Reagent: The **BCL2 (18q21) Proximal - XL** FISH probe is a direct-labeled DNA probe provided at 10 times the concentration than the amount of DNA probe used in the QC test procedure.

Kreatech FISH probes are REPEAT-FREE™ and therefore do not contain Cot-1 DNA Hybridization efficiency is increased and background, due to unspecific binding, is highly reduced.

Patterns: Two green (2G) signals will identify the non-aberrant BCL2 loci.

	Signal Pattern
Expected Signals	2G

Warning and precautions: In case of emergencies check SDS sheets for medical advice. SDS sheets may be obtained by either contacting Leica Technical Support or visiting www.LeicaBiosystems.com. DNA probes contain formamide which is a teratogen; do not inhale or allow skin contact. Wear gloves and a lab coat when handling DNA probes. All materials should be disposed of according to your institution's guidelines for hospital waste disposal.

Reagent Storage and Handling:

Store at 2-8 °C. Reagents should not be used after the expiration date on the vial label.

TECHNICAL SUPPORT

Technical support is available at www.LeicaBiosystems.com/service-support/technical-support/ or toll free at 800-248-0123 or via e-mail: kreatech-support@leicabiosystems.com.

CUSTOMER SERVICE

Kreatech probes may be ordered through Leica Customer Service toll free at 800-248-0123 or order via e-mail: purchase.orders@leica-microsystems.com.



Manufactured by
Kreatech Biotechnology B.V.
An entity of Leica Biosystems