

Kreatech™ FISH probes

Product Information Sheet

KI-10102


DLEU1 (13q14) / 13qter

100 µl

DANGER



FORMAMIDE

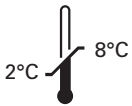


Kreatech Biotechnology B.V.
Vlierweg 20
1032 LG Amsterdam
The Netherlands
www.LeicaBiosystems.com

RUO - Research Use Only

Not for use in diagnostic procedures

PI-KI-10102_D2.2

Published August 2015

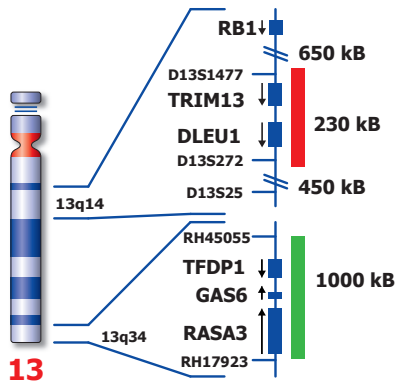

Not to scale

KI-10102

Kreatech™ DLEU1 (13q14) / 13qter FISH probe

Introduction: The **DLEU1 (13q14)** specific FISH probe is optimized to detect the DLEU1 gene region at 13q14. The **13qter (13q34)** FISH probe is included to facilitate chromosome identification.

Critical region 1(red): The **DLEU1 (13q14)** specific FISH probe is direct-labeled with PlatinumBright™550.
Control region 2(green): The **13qter** specific FISH probe is direct-labeled with PlatinumBright™495.

Reagent: Kreatech probes are direct-labeled DNA probes provided in a ready-to-use format. Apply 10 µl of probe to a sample area of approximately 22 x 22 mm.

Please refer to the Instructions for Use for the entire Kreatech FISH protocol.

Kreatech FISH probes are REPEAT-FREE™ and therefore do not contain Cot-1 DNA. Hybridization efficiency is increased and background, due to unspecific binding, is highly reduced.

Patterns: The **DLEU1 (13q14) / 13qter** FISH probe is designed as a dual-color assay to detect deletions at 13q14. Deletions involving the DLEU1 (13q14) gene region will show one red signal and two green signals at the 13qter control region (1R2G). Two single color red (R) and green (G) signals will identify the normal chromosomes 13 (2R2G).

	Normal Signal Pattern	Del(13q14)	Del (13q14-13q34)
Expected Signals	2R2G	1R2G	1R1G

References: Stilgenbauer S et al., 1998, Oncogene, 16; 1891 – 1897
Wolf S et al, 2001, Hum. Molec. Genet., 10; 1275-1285
Elinaei M et al., 2003, Genes Chromosomes Cancer, 36; 99 – 106

Warning and precautions: In case of emergencies check SDS sheets for medical advice. SDS sheets may be obtained by either contacting Leica Technical Support or visiting www.LeicaBiosystems.com. DNA probes contain formamide which is a teratogen; do not inhale or allow skin contact. Wear gloves and a lab coat when handling DNA probes. All materials should be disposed of according to your institution's guidelines for hospital waste disposal.

Reagent Storage and Handling: Store at 2-8 °C. Reagents should not be used after the expiration date on the vial label.

TECHNICAL SUPPORT Technical support is available at www.LeicaBiosystems.com/service-support/technical-support/ or toll free at 800-248-0123 or via e-mail: kreatech-support@leicabiosystems.com.

CUSTOMER SERVICE Kreatech probes may be ordered through Leica Customer Service toll free at 800-248-0123 or order via e-mail: purchase.orders@leica-microsystems.com.